

Highspeed Handpieces



Discover the advantages...



Discover the advantages





Good reasons

Good arguments for MK-dent Turbines



Robust 25.000 LUX fibre optic glass rod for perfect view in every situation. Sterilization resistant at up to 135 °C and thermally disinfectable without reducing the light intensity



Perfect working and easy to care for 3-port- or 1-port spray cooling-system. A good prevention for thermal damage of the tooth



3 different head types and sizes for perfect working in all areas. Choose between Power-Head. Small-Head and Mini-Head



Ergonomic shape with a special titanium plated casing. Easy to clean surface, nearly as hard as a diamond. "Titanium Grip" – the ingenious extra touch for the daily work



High-Tech ceramic-bearings for a quiet and vibration-free running with a long lifetime



Comfort pushbutton chuck system for quick, easy and safety bur changing



The concept of compatibility - one more good reason to decide for MK-dent



^{**} Further information's to our couplings you'll find on page 9



TITANIUM LINE (with light)





TITANIUM LINE (without light)





CLASSIC LINE (with light)





CLASSIC LINE (without light)





BASIC LINE



HC2011 - Power head, 1-port spray, ceramic-bearings, comfort-pushbutton-system, 18 Watt power at 340.000 rpm, fits for 2-hole Borden tubings



HC2012 - Power head, 1-port spray, ceramic-bearings, comfort-pushbutton-system, 18 Watt power at 340.000 rpm, fits for 4-hole tubings



HC2013 - Small head, 1-port spray, ceramic-bearings, comfort-pushbutton-system, 15 Watt power at 420.000 rpm, fits for 2-hole tubings



HC2014 - Small head, 1-port spray, ceramic-bearings, comfort-pushbutton-system, 15 Watt power at 420.000 rpm, fits for 4-hole tubings



HC1011 - Power head, 1-port spray, ceramic-bearings, friction-grip-system, 18 Watt power at 340.000 rpm, fits for 2-hole tubings



HC1012 - Power head, 1-port spray, ceramic-bearings, friction-grip-system, 18 Watt power at 340.000 rpm, fits for 4-hole tubings



Coupler





ALL-SYNTHETIC OIL

All-Synthetic Special Oil for Lubrication & Cleaning of Handpieces For any well-known Brand of Handpiece

Use for:

- Turbines
- Straight and contra angle Handpieces
- Air Motors

- Surgical Instruments
- Air Scalers

Premium Quality

MK-dent Premium Handpiece Lubricant & Cleaner is an all-synthetic special oil with the highest purity class in combination with best quality additives formulated to give your hardworking instruments a much longer lifespan

Advantages compared to standard mineral oil lubrication:

- Lower working temperature and more power because off less friction
- Less deposits
- Less consumption, more economical
- Better flow and allocation behavior
- Better wear protection, thus a longer lifespan
- Best thermal stability and oxidation resistance
- Excellent material compatibility
- Higher performance

High Performance

Keep in mind that your instruments spin up to 420.000 rpm. The parts spinning such high speed revolution need a certain oil. Protect your instruments from premature wear and tear with MK-dent Premium Lubricant and Cleaner.

Premium Handpiece Lubricant and Cleaner The ALL-SYNTHETIC special oil for any brand of handpiece CAUTION: CONTENTS UNDER PRESSURE. MAD BACK PAMEL CAREFULLY. Volume - 500ml 650 (6 NET 12.20 oz (3459)

Lubrication Tools:






Spare Parts

Spray- and Backcap keys, stainless steel



Spray- and Backcap keys for MK-dent Turbines and most of KaVo* Turbines

Gaskets and O-Rings

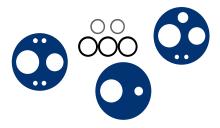

Gaskets and O-Rings for cartridges, turbine housings, spray cover and coupler for MKdent and most KaVo* turbines and couplers

Replacement Cartridges

Fits for MK-dent Turbines (incl. 2 O-Rings)					
TC1011	FG-Cartridge fits for HC1011 and HC1012 (Ceramic Bearings)	-			
TC2031	Pushbutton - Cartridge fits for Powerhead (Ceramic Bearings)				
TC2032	Pushbutton - Cartridge fits for Small- and Mini-Head (Ceramic Bearings)				

Fits for original KaVo* Turbines (incl. 2 O-Rings)					
TC1011K	FG-Cartridge fits for Super Torque 625/630 (Ceramic Bearings)				
TC2011K	Pushbutton - Cartridge fits for Super Torque 630B/C und 640B/C (Ceramic Bearings)	-0/1			
TC2650BK	Pushbutton - Cartridge fits for Super Torque 650B (Ceramic Bearings)	-100			
TC2655BK	Pushbutton - Cartridge fits for Super Torque 655B/C (Ceramic Bearings)	-070			
TS2647BK	Pushbutton - Cartridge fits for Super Torque 647B und 649B (Steel Bearings)	400			
TC2700BK	Pushbutton - Cartridge fits for Gentle Force 7000 (Ceramic Bearings)	4			
TC2012K	Pushbutton - Cartridge fits for Bella Torque 639B/C und 642B/C (Ceramic Bearings)				
TS2635BK	Pushbutton - Cartridge fits for Mira 635B (Steel Bearings)				
TS2637BK	Pushbutton - Cartridge fits for Bella Torque 637B (Steel Bearings)				

Bulbs



You can find more spare parts in our spare part catalogue, at our website **www.mk-dent.com** or contact your authorized MK-dent dealer





MK-dent is committed to sell high quality dental equipment and supplies through authorized dealer network 's worldwide. More about us and our program, including Low Speed Handpieces, Rotary Instruments, spare parts, tools, care products and more, you'll find at our website: www.mk-dent.com, in our special catalogues, data sheets and of course by personal information through our dealers and German headquarter.

Authorized Dealer:		